
OPTICELL SOFT Plus

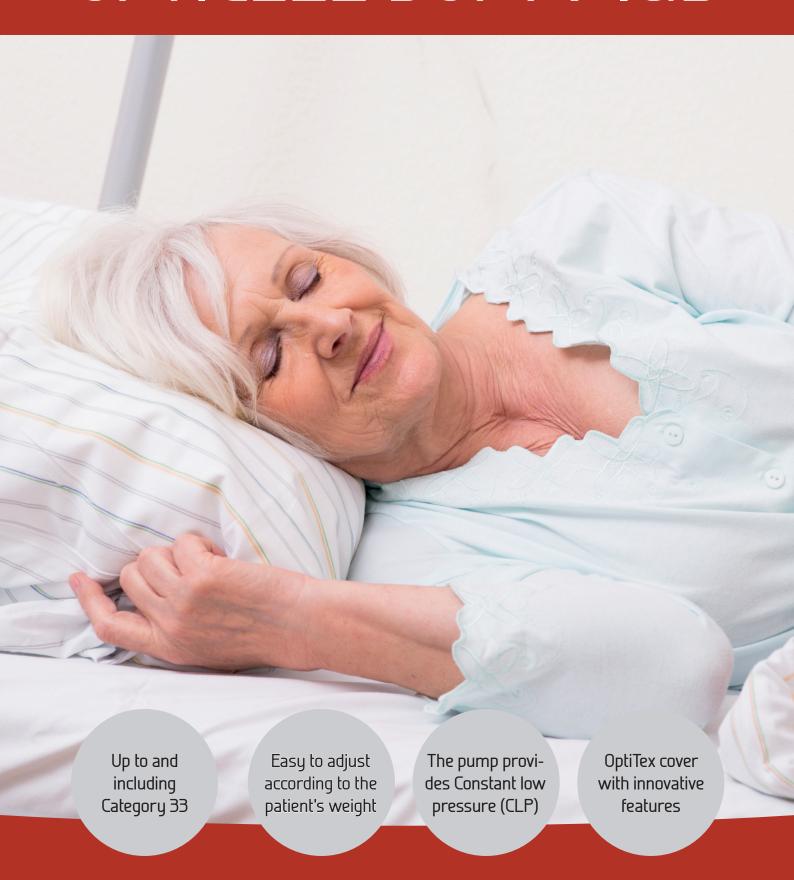

A static pressure-relieving air mattress with hygiene cover, manual pump and interchangeable cells

OPTICELL SOFT PLUS

OptiCell SOFT PLUS, is a static pressure-relieving air mattress developed as a preventive and therapeutic aid in the treatment of pressure ulcers up to and including category 3.

OptiCell Soft Plus is a bed mattress type and is placed on top of a foam mattress (e.g. one of Järven's standard mattresses). OptiCell Soft Plus comes with a durable OptiTex cover and CPR valve for emergency situations. The 13 cm thick mattress has a cover that is impermeable to liquids. OptiCell Soft Plus comes with an energy-efficient, low-noise pump, OC-15, which is easy to use. Filling the mattress normally takes about 15 minutes.

FEATURES - Pump and mattress cells

Adaptation to patient weight, 30-250 kg.

Nursing mode with automatic return to previous mode after 30 minutes. CPR valve - In emergency mode where CPR needs to be started, the CPR

valve opens to provide a quick air discharge.

Alarm function - The air pump has an alarm function with both light and

sound alarms that warn of low pressure in the mattress.

Heel relief - The mattress is equipped with a heel relief system that allows easy disconnection of the cells using an integrated air lock.

FEATURES - OPTITEX COVER

Integrated cable holder - The OptiTex cover has an integrated, liquid-tight cable holder, which prevents the spread of bacteria, keeps the power cord out of the way and protects the air hoses.

Hygienic - The OptiTex cover can be sanitised with disinfectant and can also be machine washed at 95 $^{\circ}\text{C}.$

Cells and hoses are easily removed from the cover and rinsed clean with warm water (max 70°C) or disinfectant.

The OptiTex cover allows release of air from of the cells while remaining liquid-tight. This reduces the risk of shear and other complications for the user. Drip-proof waterproof zip on four sides. All seams are welded. Underside with four elasticated straps.

Produktspecifications

STANDARD SIZES 190/200/210/220x80/85/90/105/120x13 cm. For other sizes please contact Järven Health Care.

PRESSURE ULCER CATEGORY 1-3*

CELLS Replaceable cells made of soft nylon/TPU (polyurethane).

COVER OptiTex bi-elastic double-coated textile with soft elastic PU.

PUMP OC-15. Sound level: 21 db. 230 V / 50 Hz. A Ready light indicates when the pump and mattress are ready to use. If the light is not on, it is an indication that something is wrong.

WARRANTY 2 years on manufacturing defects.

FIRE SAFETY STANDARDS SS 876 00 01, EN 597-1 & 2.

USER WEIGHT 30-250 kg.

RECYCLING Consult user manual

WEIGHT Mattress 6 kg (200 x 90 cm), pump 2,0 kg

LABELLING Date of manufacture, EAN, care instructions, fulfilled standards and placement in bed.



* Categories defined by the International Pressure Ulcer Classification system. Guidelines produced by the European and US National Pressure Ulcer Advisory Panels (EPUAP and NPUAP), together with the Pan Pacific Pressure Injury Alliance (PPPIA).



